

» i-Sleeve™

Sleeve

Supplier:

Lapoc Soft Solution (Imasen)

Features & Benefits:

- Can be used for a vacuum or suspension systems
- Airtight
- Lightweight thin, rubber based gel material provides a comfortable fit for patient
- Inside surface is treated to provide anti-bacterial deodorizing properties
- Washable
- Latex free
- Special tape allows for repair of small tears

Specifications:

- Material thickness: *1mm*
- Sizes: *X-Small, Small, Medium, Medium-Large, Large and X-Large*
- Colors: *Beige and Black*
- Distal circumference is measured at patella
- Proximal circumference is measured at mid thigh

Indications:

- Transtibial applications

Call Points:

- Prosthetist
- Prosthetic technician

Warranty:

- Against manufacturers defects




i-Sleeve™ » Lapoc Soft Solution (Imasen)

Item No.	Description	Size	Patella Circ.	Thigh Circ.	Color
IEC-M1601BGE-SS	Suspension sleeve	X-Small	27–30 cm (10 ⁵ / ₈ "–11 ³ / ₄ ")	38–43 cm (15"–17")	Beige
IEC-M1601BLK-SS	Suspension sleeve	X-Small	27–30 cm (10 ⁵ / ₈ "–11 ³ / ₄ ")	38–43 cm (15"–17")	Black
IEC-M1601BGE-S	Suspension sleeve	Small	30–35 cm (11 ³ / ₄ "–13 ³ / ₄ ")	43–50 cm (17"–19 ³ / ₄ ")	Beige
IEC-M1601BLK-S	Suspension sleeve	Small	30–35 cm (11 ³ / ₄ "–13 ³ / ₄ ")	43–50 cm (17"–19 ³ / ₄ ")	Black
IEC-M1601BGE-M	Suspension sleeve	Medium	34.5–37.5 cm (13 ⁵ / ₈ "–14 ³ / ₄ ")	50–56.5 cm (19 ³ / ₄ "–22 ¹ / ₄ ")	Beige
IEC-M1601BLK-M	Suspension sleeve	Medium	34.5–37.5 cm (13 ⁵ / ₈ "–14 ³ / ₄ ")	50–56.5 cm (19 ³ / ₄ "–22 ¹ / ₄ ")	Black
IEC-M1601BGE-ML	Suspension sleeve	Medium-Large	35–38 cm (13 ³ / ₄ "–15")	56–63 cm (22"–24 ⁷ / ₈ ")	Beige
IEC-M1601BLK-ML	Suspension sleeve	Medium-Large	35–38 cm (13 ³ / ₄ "–15")	56–63 cm (22"–24 ⁷ / ₈ ")	Black
IEC-M1601BGE-L	Suspension sleeve	Large	36–42 cm (14 ¹ / ₄ "–16 ¹ / ₂ ")	56–63 cm (22"–24 ⁷ / ₈ ")	Beige
IEC-M1601BLK-L	Suspension sleeve	Large	36–42 cm (14 ¹ / ₄ "–16 ¹ / ₂ ")	56–63 cm (22"–24 ⁷ / ₈ ")	Black
IEC-M1601BGE-XL	Suspension sleeve	X-Large	37–44 cm (14 ¹ / ₂ "–17 ³ / ₈ ")	61–68 cm (24"–26 ³ / ₄ ")	Beige
IEC-M1601BLK-XL	Suspension sleeve	X-Large	37–44 cm (14 ¹ / ₂ "–17 ³ / ₈ ")	61–68 cm (24"–26 ³ / ₄ ")	Black

i-Sleeve™ » Repair Information

STEP.1


Set the temperature of iron
150 – 160 degrees



※ Note: Do not use steam


STEP.2

Cut repair tape little larger than the size of the hole, with round edges

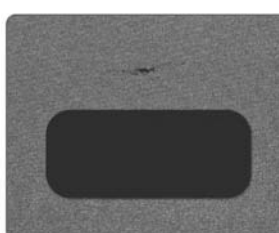


STEP.3

Set i-Sleeve on a flat repair surface

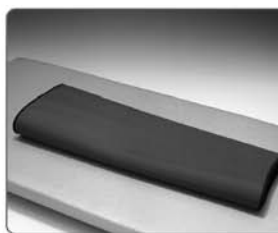


Set the repair tape on the hole




STEP.4

Set the wax paper on i-Sleeve



STEP.5


When at temperature, press iron onto the wax paper, over the tape, for 20 seconds



※ Note: Be careful not to burn

STEP.6

Make sure repair tape attached to the i-Sleeve



If repair tape does not attach to the i-Sleeve, repeat **STEP 3** again.

